Anabolic Lab

Analytical Lab Testing for Anabolic Steroid Harm Reduction

https://anaboliclab.com

Lab Report – May 27, 2016

Pharmacom Labs PHARMA Tren H100 (100mg/mL) Trenbolone Hexahydrobenzylcarbonate



If you like the attached report, please visit <u>our website</u> to find out how you can support independent and objective lab testing.



Simec AG

Areal Bleiche West Postfach 413 CH - 4800 Zofingen Telefon +41 62 752 83 08 Fax +41 62 752 83 09 E-Mail: info@simec.ch www.simec.ch

Test report

Number

C 42369

Object of

analysis:

Pharmacom Labs Pharma Tren H100

Customer:

AnabolicLab.com

Samples: 1

Batch:

604331, 25a7j3exuu

Receipt: 24.03.2016

Subject:

Assay, Identification Trenbolone hexahydrobenzylcarbonate

Method(s):

HPLC

Sample / Analysis:	Method	Status	Result	Specification	complies
Trenbolone hexahydro- benzylcarbonate: Identification Content	HPLC HPLC	N N	complies 94.8 mg/mL	complies 100 mg/mL	yes -
TAMC TYMC		N N	< 1 cfu/mL < 1 cfu/mL	*	-

Method-status:

G = GMP A = accreditated V = generally validated P = validated on product N = not validated

E = external lab

Remarks:

The sample(s) mentioned above have been analysed as they have been sent to us by the client. We have not controlled external processes like manufacturing, labelling, sampling, shipping, storage. These results are for information only and do not compensate for correct quality control by the manufacturer/distributor. Generally, pharmaceutical products have to be produced and distributed under full GMP/GDP regime, including GMP-compliant analyses with validated

methods.

Signatures:

controlled & released:

Staff:

completed: 27.05.2016

The enclosed results refer exclusively to the object of examination described above. The measuring accuracy is available on request. Only handwritten signatures are valid. Without written permission, it is not allowed to publish single parts of this report.

page 1 of 1

Simec AG Page 1 of 1 (38)

Druckdatum:

Injektvolumen:

17.05.2016 08:35:10

20

Report steroid analysis

HPLC:

Nr.5

Analysendatum:

14.05.2016 18:59:02

Vial:

Probenidentifikation:

167 C42369 A

Operator:

admin

Chrom.-Beding.

ACN / Wasser, gradient

Temperatur:


40 °C

Säule

Nr. 121, 250 x 4,6, C18

Detektion:

UV 240nm

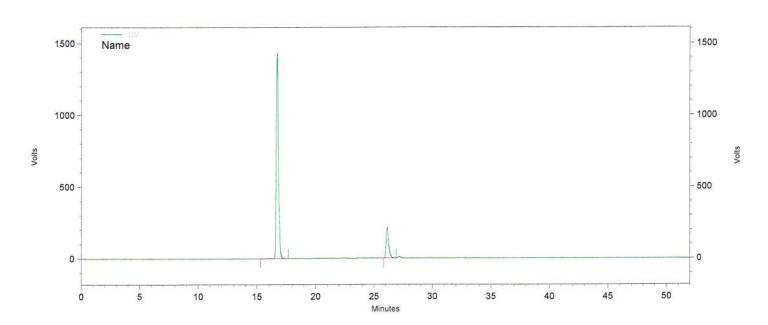

Channel A

UV	Resu	ts
----	------	----

Name	Time	Area	Amount	Area Percent
Trenbolone acetate	2000000		0,000 BDL	
	16,707	73595879	0,000	85,174
Nandrolone Phenylpropionate	8		0,000 BDL	
,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,	26,117	12810408	0,000	14,826
Trenbolone			0,000 BDL	
Hexahydrobenzylcarbonat			72	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
Totals				
		86406287	0,000	100,000

Daten File:

C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse

Steroide\125_C42369 A-Rep1.dat

Methode:

Simec AG Page 1 of 1 (39)

Report steroid analysis

HPLC: Nr.5

Analysendatum: 14.05.2016 19:52:21 Druckdatum: 17.05.2016 08:35:14

Vial: 167 Injektvolumen: 20


Probenidentifikation: C42369 A Operator: admin

Chrom.-Beding. ACN / Wasser, gradient

Temperatur: 40 °C

Säule Nr. 121, 250 x 4,6, C18

Detektion: UV 240nm

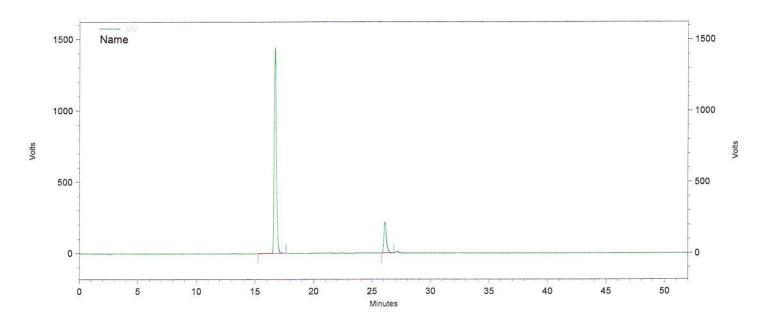

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Trenbolone acetate			0,000 BDL	
	16,687	73763959	0,000	85,145
Nandrolone Phenylpropionate	**************************************		0,000 BDL	
	26,087	12869687	0,000	14,855
Trenbolone			0,000 BDL	
Hexahydrobenzylcarbonat			18°	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
Totals				100 January - Mars 2000
		86633646	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse

Steroide\125_C42369 A-Rep2.dat

Simec AG Page 1 of 1 (40)

Report steroid analysis

HPLC: Nr.5

Analysendatum: 14.05.2016 20:45:42 Druckdatum: 17.05.2016 08:35:19
Vial: 167 Injektvolumen: 20


Vial: 167 Probenidentifikation: C42369 A

Operator: admin
Chrom.-Beding. ACN / Wasser, gradient

Chrom.-Beding. ACN / Wass Temperatur: 40 °C

Säule Nr. 121, 250 x 4,6, C18

Detektion: UV 240nm

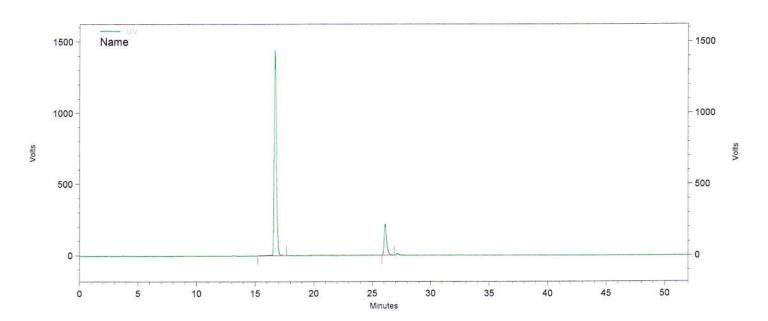

Channel A

UV Results Name	Time	Area	Amount	Area Percent
Trenbolone acetate			0,000 BDL	
	16,677	73888528	0,000	85,149
Nandrolone Phenylpropionate			0,000 BDL	The Control of the Co
Harlarolone Frienty propionate	26,097	12886543	0,000	14,851
Trenbolone			0,000 BDL	9
Hexahydrobenzylcarbonat			9.50)	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
			2004-00-00-00-00-00-00-00-00-00-00-00-00-	
Totals				
		86775071	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse

Steroide\125_C42369 A-Rep3.dat

Simec AG Page 1 of 1 (41)

Report steroid analysis

HPLC: Nr.5

Analysendatum: 14.05.2016 21:39:08 Druckdatum: 17.05.2016 08:35:23
Vial: 168 Injektvolumen: 20

Vial: 168
Probenidentifikation: C42369 B


Operator: admin

Chrom.-Beding. ACN / Wasser, gradient

Temperatur: 40 °C

Säule Nr. 121, 250 x 4,6, C18

Detektion: UV 240nm

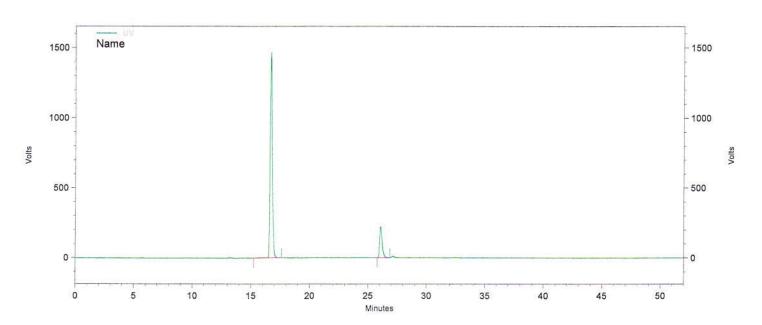

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Trenbolone acetate			0,000 BDL	
	16,677	74682628	0,000	85,078
Nandrolone Phenylpropionate	\$100 P. \$100 S. \$100 S		0,000 BDL	*** ##********************************
	26,080	13098507	0,000	14,922
Trenbolone	50		0,000 BDL	52
Hexahydrobenzylcarbonat			·00	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
Totals				
		87781135	0,000	100,000

Steroide\126_C42369 B-Rep1.dat

Simec AG Page 1 of 1 (42)

Report steroid analysis

HPLC: Nr.5

Analysendatum: 14.05.2016 22:32:24 Druckdatum: 17.05.2016 08:35:28 Vial: Injektvolumen: 20

Vial: 168 Injektvolumen: Probenidentifikation: C42369 B


Operator: C42369 B

Chrom.-Beding. ACN / Wasser, gradient

Temperatur: 40 °C

Säule Nr. 121, 250 x 4,6, C18

Detektion: UV 240nm

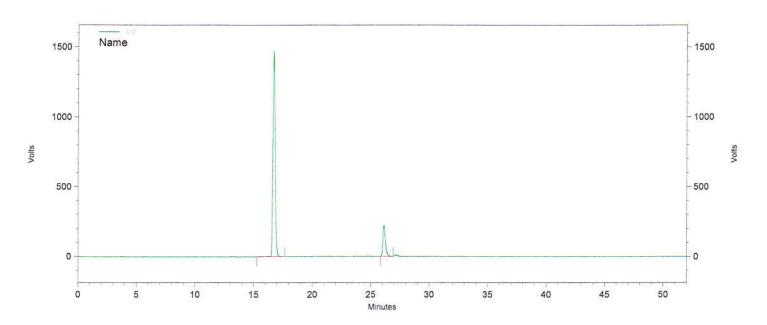

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Trenbolone acetate			0,000 BDL	
	16,680	75309950	0,000	85,081
Nandrolone Phenylpropionate	***************************************		0,000 BDL	
20 24 1958	26,117	13205676	0,000	14,919
Trenbolone			0,000 BDL	19
Hexahydrobenzylcarbonat			ů.	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
Totals				
		88515626	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse

Steroide\126_C42369 B-Rep2.dat

Simec AG Page 1 of 1 (43)

Druckdatum:

Injektvolumen:

17.05.2016 08:35:32

20

Report steroid analysis

HPLC:

Nr.5

Analysendatum:

14.05.2016 23:25:49

Vial:

168

Probenidentifikation:

C42369 B

Operator:

admin

Chrom.-Beding.

ACN / Wasser, gradient

Temperatur:


40 °C

Säule

Nr. 121, 250 x 4,6, C18

Detektion:

UV 240nm

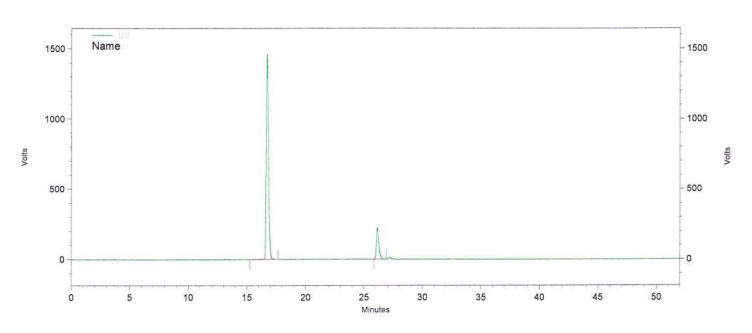

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Trenbolone acetate			0,000 BDL	
	16,697	75668869	0,000	85,128
Nandrolone Phenylpropionate	0000 € 0000000		0,000 BDL	
	26,133	13219689	0,000	14,872
Trenbolone	C1990 900 m ■ 000 0399 (00 000 00		0,000 BDL	40 600m25 #34 0 40094 1 0 m25
Hexahydrobenzylcarbonat			**************************************	
Testosterone enanthate			0,000 BDL	
Testosterone cypionate			0,000 BDL	
Totals				
		88888558	0,000	100,000

Daten File:

C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse

Steroide\126_C42369 B-Rep3.dat

Methode: